


Reach the next level of LabX titration

Version 1.x upgrade

Series	State
pH Electrode ...	done
pH Electrode ...	done
pH Electrode ...	done
pH Electrode ...	done
pH Electrode ...	done
KHP Content	active
KHP Content	ready

Method ID	Title
M560	pH-Determination in W
M561	Acid Content in Wines
M562	Combined pH and Acic
M563	Free Sulphur Dioxide i
M564	Total Sulphur Dioxide
M566	Determination of the f
M567	Determination of the f
M568	Reducing Sugar accor
M569	Ascorbic Acid (Vitamin
M570	Ascorbic Acid (Vitamin
M573	Ascorbic Acid (Vitamin

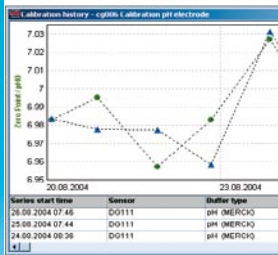

Simpler operation

In addition to generally improved user interface functionality such as context sensitive right-click menus and mouse-over help texts, LabX titration offers simplified management of scheduled titrations: The new job queue displays all started methods and series and allows the user to edit, stop and delete these processes.

Greater flexibility

Big improvements that make a big difference:


- Integrated application database with 500+ methods
- Individual report layout
- User defined data organization
- Predefineable filters
- More powerful LIMS integration

Highest security

The high security level is further enhanced to achieve full support of compliance with FDA 21 CFR part 11: Extended user management functionality, total instrument management, graphical filter and calibration history as well as expanded validation documentation are now available.

More power with Titration Excellence

In connection with Excellence Titrators, LabX titration offers many additional features such as Dual Mode (work at the PC and at the titrator simultaneously), One Click™ Titration (titration starts with just one click) and a much improved method and series data editor.

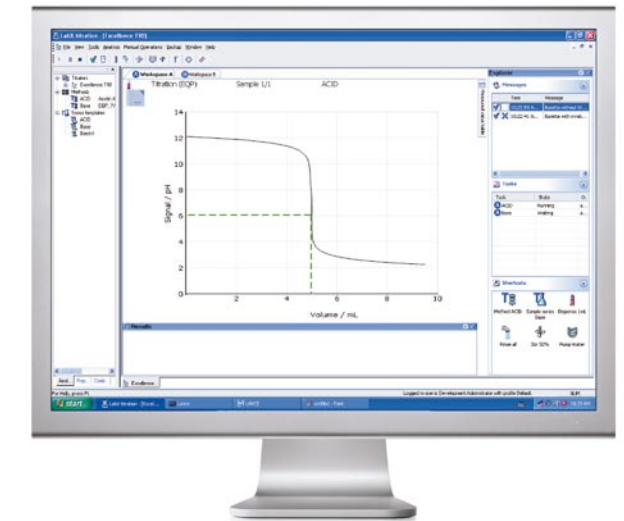

Gain flexibility and power with LabX titration 2.6

LabX titration PC software offers data security, traceability, ease of use and efficient operation. Version 2.6 now includes many improvements based on the feedback of thousands of users. It fully supports the current DL3x, DL5x and DL7x titrators as well as the new Titration Excellence line.

Main new features of LabX titration version 2.6 compared to 1.x

User Interface improvements		pro	light
Integrated application database with 500+ methods	Directly search for an appropriate method for your application within LabX, read the application description and use the method right away.	•	•
User defined data organization	Create your own structure for reports storage. However, the elegant automatic structure is still available anytime – just in case you miss something.	•	
Add notes to each sample	Add your notes to each sample. The notes will be available in reports as well as on printouts and in exported data sets.	•	•
Predefineable filters	Save the filters you often use, e.g. your own results since the beginning of the month.	•	•
New languages	English, German, French and now Chinese, Spanish, Italian and Russian.	•	•
Expanded mouse actions	Improved right click functionality as well as mouse-over-information and tool tips.	•	•
Greater flexibility		pro	light
Use very long and multiple sample IDs	LabX allows sample IDs with a total length of up to 120 characters or multiple IDs with a total length of 120 characters.	•	•
Job queue	You may start as many series of single measurements as you want. The job queue window will always show you what's done and what's coming next and is fully editable. Change the order of tasks in the job queue or edit items in the job queue.	•	•
User defined data	Import/input results from any other source such as other instruments or from manual measurements and use them to perform additional calculations.	•	
Individual reports	Define your reports with the layout editor, e.g. add your company logo, customize headers and footers etc.	•	
Export and import of report templates	Your report templates can now be imported and exported.	•	•
Direct pdf generation	Create quickly a pdf of any report in LabX.	•	•
More flexible file export	Any delimiter, file extension etc. can be directly set within LabX. Therefore LabX can now easily generate csv and custom file formats.	•	
No licensing fees for additional PCs	The new licensing system of LabX pro titration is based on the number of titrators you use with LabX. This means you can now install LabX pro on as many computers as you want without increasing your license fee.	•	
Directly connect titrators over Ethernet	You do not need to put a PC close to each titrator. There is no serial connection to the PC necessary anymore!	•	•
Higher security		pro	light
Extended security and user management functionality	Additional policies for password protection and an even more powerful user manager make LabX more secure than ever.	•	
Access violation alert	The administrator receives a message if a user account has been locked e.g. if a logon failed three times.	•	
Graphical titer and calibration history	Keep track of your calibration and titer determinations with automatic control charts.	•	•
Expanded validation documentation	The LabX validation handbooks I & II have been expanded.	•	



Quality certificate ISO 9001
 Environment certificate ISO 14001
 Internet: <http://www.mt.com>
 Worldwide service

Subject to technical changes
 ©08/2007 Mettler-Toledo AG
 Printed in Switzerland 51724594

www.mt.com/LabX

Visit for more information